




For classic fat extraction

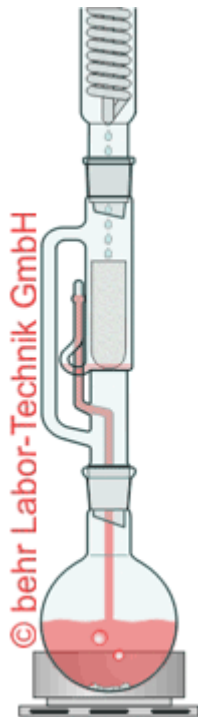

Complete compact extraction system with basic rack, heating device, holder, hoses and glassware (reaction flasks, extractor, Dimroth condenser for extraction). Infinitely variable heating regulation After the extraction cycle, the compact system removes the solvent, by means of a discharge valve on the extractor directly into the reservoir.

Includes sample package, extraction sleeves and boiling stones.

Behr equipment for Soxhlet extraction comply with the most various requirements of daily laboratory routines.

With Behr hydrolysis units, acid digestion is possible before extraction (determination of total fat content according to Weibull and Stoldt).

- Practical holders for condenser and extraction adapter for secure placement between extractions
- Extractor sizes from 30 ml to 1000 ml
- Condenser with screw connections